

Overview

This repository provides wrapper code to use the generalist cell/nuclei segmentation package [Cellpose](#).

Note, this document is also provided as [PDF](#).

Two main functionalities are provided

1. Prepare z-stacks for segmentation by performing z-stack projections.
2. Provide convenient interface to batch segment a large number of images.

Table of contents

- [Overview](#)
- [Installation](#)
 - [Python with Miniconda](#)
 - [Create dedicated environment to run Cellpose](#)
 - [Installing Cellpose](#)
 - [Installing this package](#)
- [Data](#)
 - [Data organization](#)
 - [Test data](#)
- [Usage](#)
 - [Working with Jupyter notebooks](#)
 - [What is Jupyter?](#)
 - [Starting a Jupyter notebook](#)
 - [Preprocessing](#)
 - [Performing segmentation of cells and nuclei](#)
 - [Post-processing](#)
 - [Create FISH-quant \[Matlab\] outline files](#)

Installation

The installation consists of the following steps

1. **Python**. Recommended with Miniconda:
2. Create a **dedicated environment** with Jupyter to run your code.
3. Install **CellPose**.
4. Install **code from this repository**.

Python with Miniconda

We recommend installing an [Miniconda distribution of Python](#): choose Python 3.7 and your operating system.

We then recommend using the `anaconda` prompt that is available to execute the different commands listed below. This guarantees that the necessary terminal scripts are available.

Create dedicated environment to run Cellpose

We recommend creating a dedicated environment to run Cellpose. To create an environment called **cellpose**, open an anaconda prompt and type. Note that you will also install jupyter, which will allow to run the jupyter notebooks for easier execution (confirm with **y** when asked if you want to proceed):

```
conda create --name cellpose python=3.7 jupyter
```

Then activate the **cellpose** environment (Note you will always have to run this command when using this workflow):

```
conda activate cellpose
```

Installing Cellpose

From your **cellpose** environment, install Cellpose and its dependencies with

```
pip install cellpose --upgrade
```

Installing this package

From your **cellpose** environment, install this package its dependencies with

```
pip install git+https://github.com/muellerflorian/segmentation/develop --upgrade
```

TODO: update to master branch once merged.

Data

Data organization

We strongly recommend the following data-organization on which this workflow has been tested.

1. Images are store as single-channel multi-z-stack tif files, e.g on tif per position and channel.
2. All raw 3D images are stored in a folder **acquisition**
3. All analysis results are stored in subfolder **analysis**, where each analysis step has a separate subfolder.

The organization of the provided test data is the following

```
├─ example_data/  
│   └─ acquisition # Folder with raw data
```

```

| | | | test_pos001_cy3.tif
| | | | test_pos002_dapi.tif
| | | | test_pos002_cy3.tif
| | | | test_pos002_dapi.tif
| | | | analysis # Folder with all analysis results
| | | | segmentation-input # Folder with projected images for
segmentation
| | | | | img-prop__test_pos001_cy3.json # json file with image properties
| | | | | test_pos001_cy3.png # Projected image
| | | | | ....
| | | | | segmentation-results # Folder with segmentation results
| | | | | test_pos001_cy3.tif
| | | | | test_pos001_cy3.tif
| | | | | test_pos001_cy3.tif
| | | | | ....
| | | | | ....
| | | | | FQ_outline # [Optional] FQ outlines
| | | | | test_pos001_cy3_outline.txt
| | | | | ....

```

Test data

Already processed test data can be downloaded from [Dropbox](#). With these data, the different analysis steps can be verified.

Usage

Working with Jupyter notebooks

What is Jupyter?

Jupyter notebooks provide a convenient way to execute code with minimal user input. We further provide notebooks with interactive control, to facilitate usage.

Code is divided into code-cells (two in the example below). The currently enabled cell (the first one in the example below) is shown with a blue frame. It can be executed by pressing **SHIFT+ENTER**

```

In [1]: # Load segmentation widget
        %run segmentation_widgets.ipynb

In [ ]: # Start interface
        ui_segment_cells_nuclei

```

Starting a Jupyter notebook

1. Open an anaconda prompt.
2. Activate the **cellpose** environment:

```
conda activate cellpose
```

3. Start Jupyter notebook:

```
jupyter notebook
```

This will launch a new browser window (or a new tab) showing the Notebook Dashboard, a sort of control panel that allows (among other things) to select which notebook to open. Here, you can then start a jupyter notebook (files ending with **.ipynb**)



4. Note that closing the browser (or the tab) will not close the Jupyter Notebook App. To completely shut it down you need to close the associated anaconda terminal.

Preprocessing

Segmentation is done on 2D images. In this step, 3D images are transformed into 2D images by applying a projection.

This is done with the jupyter notebook **preprocessing.ipynb**

1. Running the first code cell will import the user-interface.
2. Running the second code cell will display the user-interface.

Note that you have to perform this projection for each channel-type. This allows to use different projection methods for a channel.

Here the following parameters can be set:

Option	Type	Default	Description
Data path	str		Full path to folder containing data to be segmented.


Option	Type	Default	Description
Results path	str		Full path to folder where results should be stored.
Channel str	str	dapi	Unique string to identify channel that should be processed.
Projection	str	mean	Different projection types: <code>mean</code> , <code>max</code> , <code>indiv</code> . The option <code>Indiv</code> implies that a z-stack is split into individual slices, stored in subfolder for each image.

- Pressing the button `Pre-process data` will start the segpre-processing. Progress can be monitored in the tab `Log`.
- Once the segmentation is finished, results can be inspected in the lower part of the interface. The dropdown menus allow to inspect the results for cell and nuclear segmentation.

Pre-processing is performed ... for more details see tab 'Log'.
Processing finished.

Pre-processing is performed:

file ▼

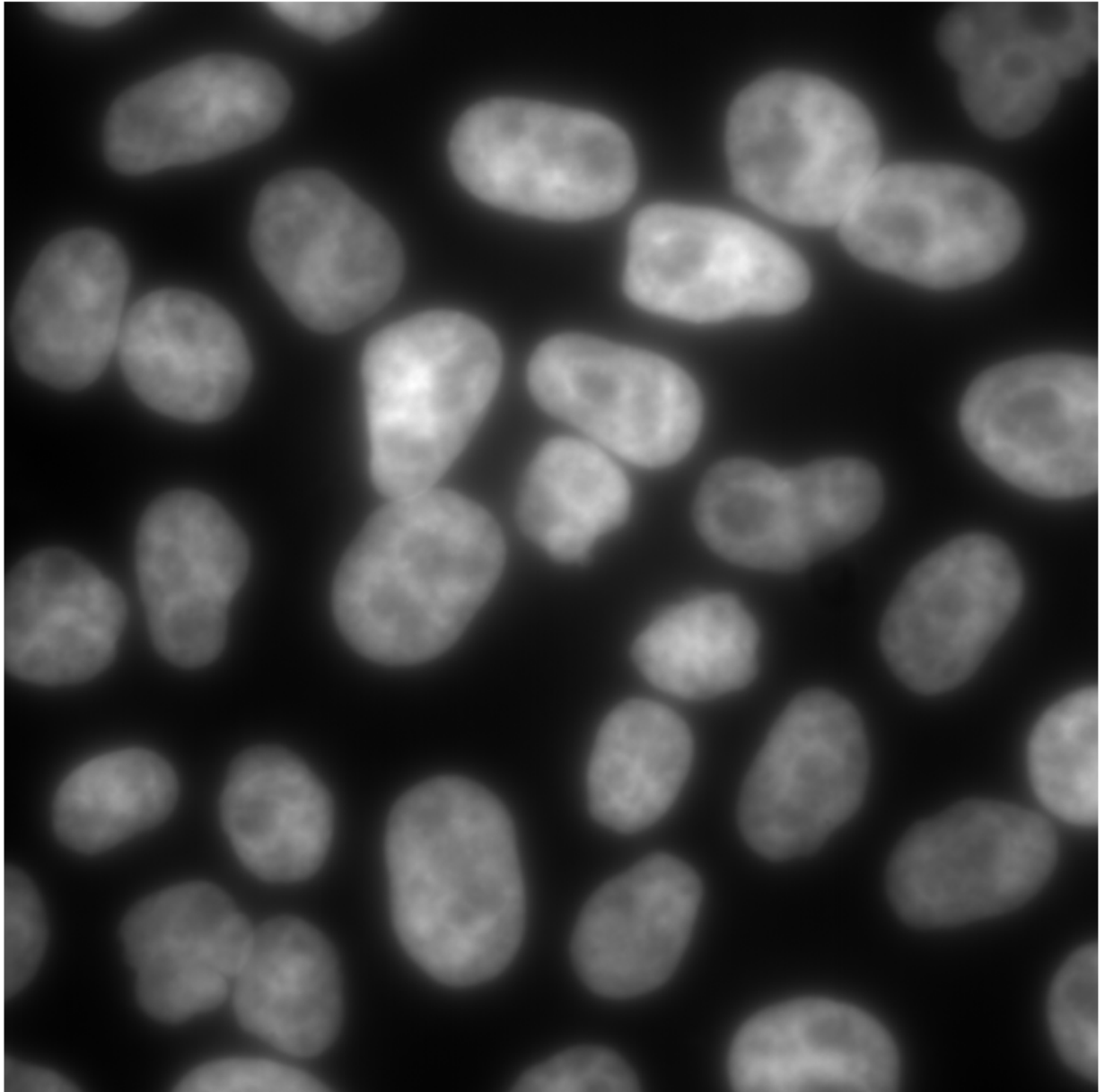

5. **Results** will be saved in the specified folder. For each image a json file with basic properties of the file, and an image with the same name as the original one will be saved.

Performing segmentation of cells and nuclei

Resizing to speed up prediction: segmentation speed depends on the image size. In our experience, resizing the images can lead to a substantial speed-up. In case you resize the images, we implemented a post-processing routine that will resize the predicted masks back to the original image size.

This is done with the jupyter notebook [segmentation_cells_nuclei.ipynb](#)


1. Running the first code cell will import the user-interface.
2. Running the second code cell will display the user-interface.

Here the following parameters can be set:

Option	Type	Default	Description
Data path	str		Full path to folder containing data to be segmented.
Results path	str		Full path to folder where results should be stored.
String cells	str	cy3	Unique identifier for images of cytoplasmic stain.
String nuclei	str	dapi	Unique identifier for images of nuclear stain.
New size	str	512, 512	String to specify new size of image. No resizing if empty.
Image ext	str	.png	File extension of images that should be segmented.
Size cells	int	100	Typical size of a cell (in resized image).
Size nuclei	int	50	Typical size of a nucleus (in resized image).

- Pressing the button **Segment data** will start the segmentation. When using CellPose for the first time, the models for nuclear and cytoplasmic segmentations are downloaded.


The actual segmentation can take a while, depending on the number of images that should be segmented (and their size). Progress can be monitored in the tab **Log**.

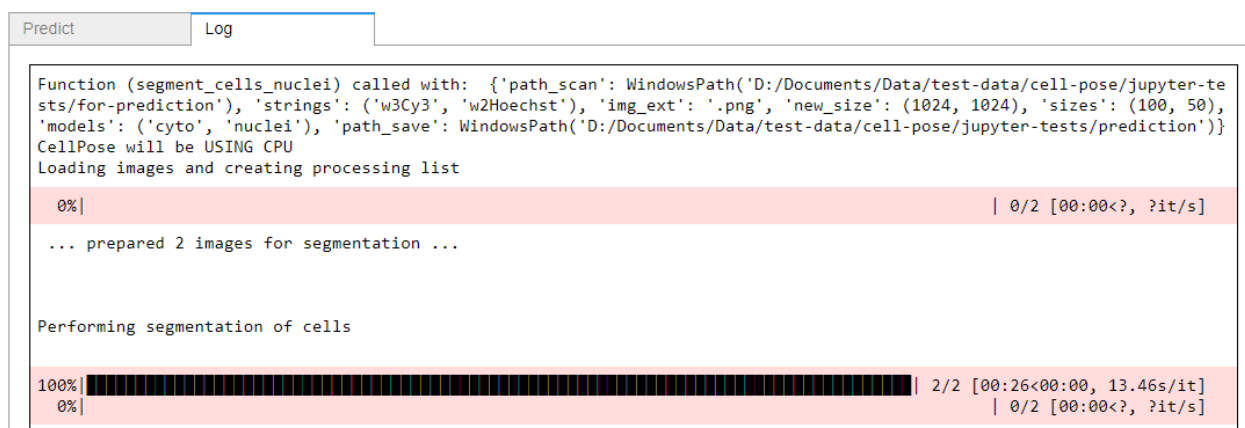


4. Once the segmentation is finished, results can be inspected in the lower part of the interface. The dropdown menus allow to inspect the results for cell and nuclear segmentation.

Performing segmentation ... for more details see tab 'Log'.
Segmentation finished

Segmentation of cells:

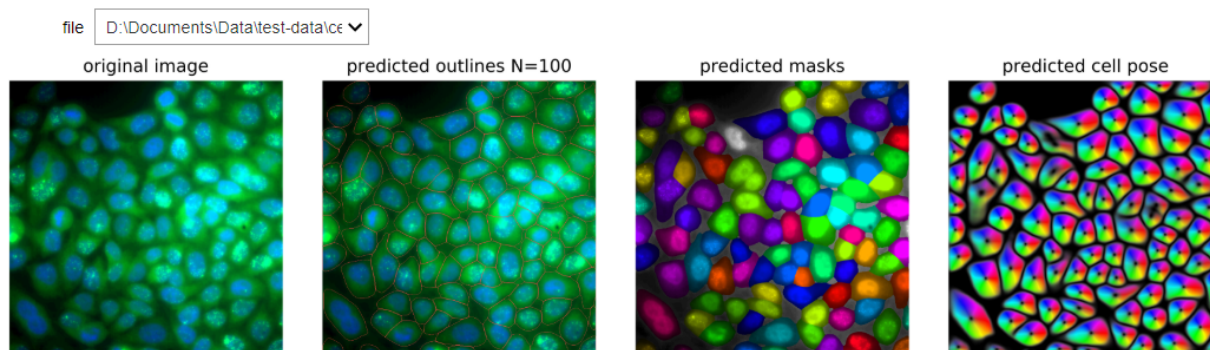

Segmentation of nuclei:

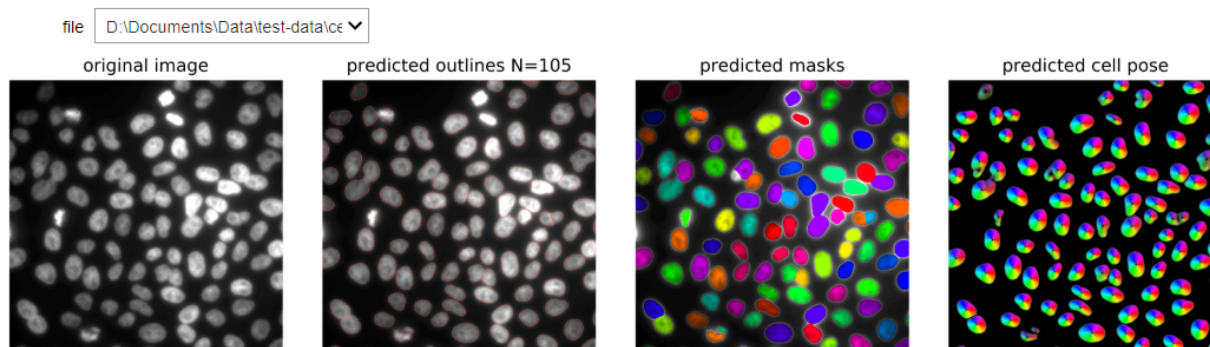

5. **Results** will be saved in the specified folder. For each image the following files, results files with different suffices are created:

- **flow_...**: these are the predicted distance maps of CellPose. They are an intermediate result, and not needed for most end-users.
- **mask_...**: these contain the actual segmentation results. Each segmented cell or nuclei is a filled object with a constant pixel value. If the images were resized during segmentation, the mask is scaled back up to the original image size. The actually obtained (smaller) mask is saved under the name **mask__rescale_...**
- **segmentation_...**: summary plot showing the input image, the predicted distance map, and the segmented objects. This plot is also shown in the interface.

Post-processing

Create FISH-quant [Matlab] outline files

In order to use the segmentation results in the Matlab version of FISH-quant, FISH-quant outline files have to be created from the mask images.

For this, we provide a Matlab GUI distributed with the [FISH-quant package](#). After installing FISH-quant, you can open this GUI from the command window with **FQ_seg**. The relevant part of the interface is the central panel *Cell Profiler: generate FQ outlines from segmentation*.

Cell profiler: generate FQ outlines from segmentation

Define experimental parameters

Identifier for images of nuclei:

Identifier for images of:

Identifier for nuclei:

Identifier for cell:

File extension of original images:

File extension of segmentation:

Part of file-name (projection) to:

☐ Generate outline files for 2nd c...

☐ Don't generate for 1st color

How to save outlines

☐ Same folder ☒ Sub-folder

☐ Generate folder

☐ Search files recur...

1. **Specify experimental parameters.** These parameters have to be defined – even if the default parameters are good. Only then the button to generate the outlines will be enabled.
2. **Define naming scheme of original images and segmentation results.** You have to define a few parameters regarding the naming convention of your files.
 1. Unique identifier for the FISH and DAPI images. These identifiers have to be defined in a way that when you take the full file-name of the FISH image and you replace the identifier for FISH, e.g. **cy3**, by the identifier of the DAPI, e.g. **dapi**, you obtain the DAPI file-name.
 2. Then, you have to define the identifier of your segmentation results. For this workflow, **mask__nuclei__** for the nuclei and **mask__cells__** for the cells.
 3. File-extensions of the masks and original file-names, e.g. **tif** for the example data.
3. [Optional] **Generating outline files for a second color.** This option allows to generate outline files for a second color, e.g. for a dual-color FISH experiment or if the cell segmentation was performed with a different channel than the FISH channel. The outlines for this color will be based on the segmentation results of the first color and the exact cells will be used. This allows a simple comparison between the detection results. As above, the identifier for the second color has to be specified, e.g. **cy5**. You also have to redefine the experimental parameters (most often to adjust the excitation and emission wavelength).

Additionally, you can choose to not create the outlines for the channel that was used for cell segmentation. This option is useful if the first color does not contain actual smFISH data but results of a dedicated cell segmentation stain.

4. Several options exist to **specify the folder** where the results will be saved. By default, the FQ outlines will be stored in a sub-folder `__FQ_outlines` within the folder containing the segmentation results. You can then move it to another location, e.g. directly into the `analysis` folder as done for the example data.
5. Specify **images that will be analyzed**. You can either choose different images that you want to analyze (`Define images`), or select an entire folder (`Select folder`). For the latter, you can also specify a recursive search; this means that all subfolders will be searched as well (not recommended for this workflow). The script will only consider images that follow the above explained naming convention – other images will be ignored. In the example data, the folder `analysis\segmentation-results` contains the relevant data.
6. **Create outlines**. Lastly, press the button `Create FQ outlines`. The script will then automatically search for the files describing the segmentation of cells and nuclei. For each image an outline file with the reference to the ORIGINAL 3D image will be generated and nuclei assigned to their respective cells.

Licensing

- **Cellpose** is (general) are BSD-licenced (3 clause)

Copyright © 2020 Howard Hughes Medical Institute

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution. Neither the name of HHMI nor the names of its contributors may be used to endorse or promote products derived from this software without specific prior written permission. THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT OWNER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

- The rest of the code provided in this repository (general) are BSD-licenced (3 clause):

Copyright © 2020, Florian Mueller
All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

- Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
- Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
- Neither the name of the copyright holder nor the names of its contributors may be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL COPYRIGHT OWNER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.